



Creative Item For Hospital and Pharmacy Innovation


CRETEM

MTC-30 (Multi Tablet Counter) Vial Packager




1. CRETEM MTC-30 Product Information

Product Name	CRETEM MTC-30 (Multi Tablet Counter) Product
Cassette Quantity	30 Cassettes
Product Structure	Standard (30 Cassettes) & Extendable to max 90 Cassettes (2 more Extender Units)
Product Control	Microsoft Surface Go or Touch Panel (Window 10) → Optional Device
Communication Method	RS422
Medication Detection	Infrared Photo-detection
Medication Storage	Cabinet type with cassettes
Dispensing Speed	Maximum Dispensing Speed: 290 pills / minute
Packing	Direct drop into cup or user's bottle
Dimension (W * D * H)	808mm * 180 mm * 1,220 mm (31.8 inch * 7.1 inch * 48 inch) / unit
Weight	60 Kg (132 lb)
Power Consumption	100W (Max. 500W)
Power Supply	220V / 60 Hz (Korea) / 120V / 60 Hz (USA) / 230V-240V / 50 Hz (Europe)
CONTACT No.	Tel No. : +82-70-7123-8522, Mobile No. : +82-10-3452-0312 (WhatsApp Available) E-mail: johnkim@cretem.co.kr

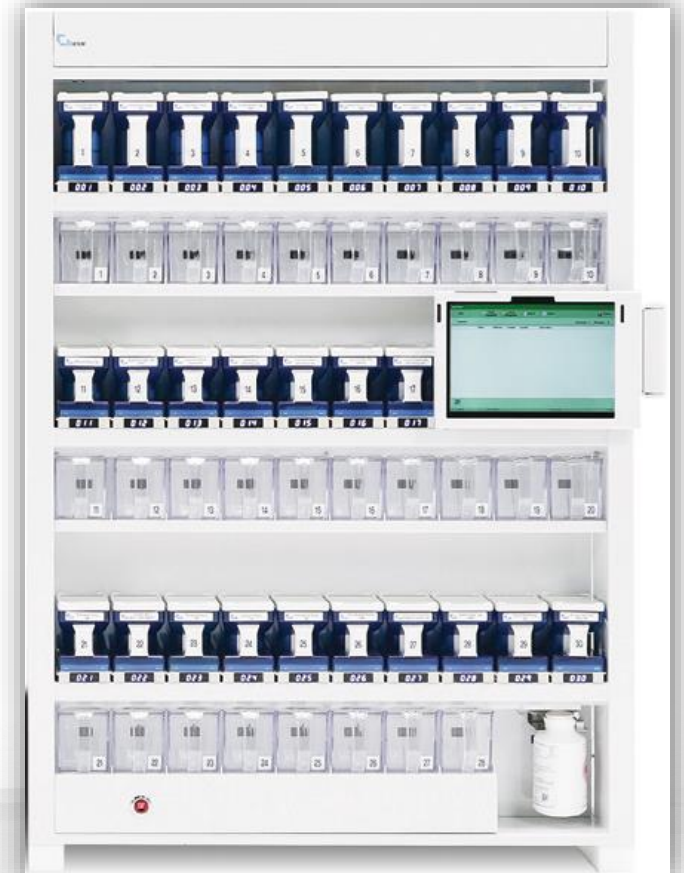

2. CRETEM MTC-30 Product Components

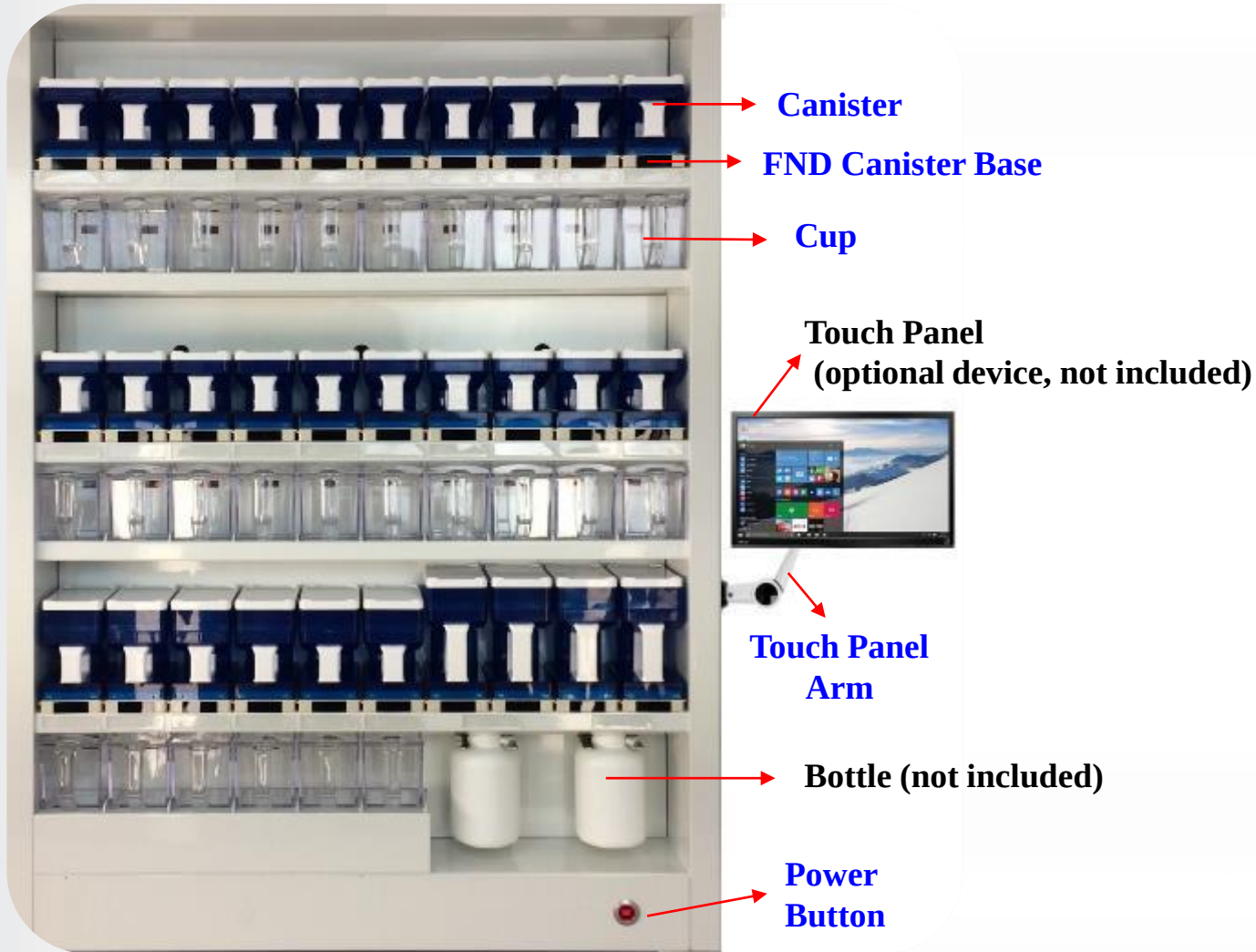

1. It can dispense medication into a bottle or cup by receiving prescription (order) in system simultaneously to enhance the pharmacy efficiency.

※ OCS (Order Communication System) or HIS(Hospital Information System) or PMS(Pharmacy Management System)



AP Counter (Operating SW on Touch Panel) ↔ MTC-30

2. Key Features

- 1) Multi-Counter for Bottle Packing of medications at same time
- 2) Control of 30 Canisters, Extendable to Max. 90 Canisters.
- 3) Compact Size for Table-Top
- 4) Fast-Dispensing within 30 seconds

3. CRETEM MTC-30 Operation Process

① Prescription (order) received from OCS / HIS / PMS.



② Put the cup into the position / Send the order to MTC.



③ MTC starts multiple orders at the same time. ("Dispensing" displayed on touch panel)



④ Printed out labels at the order in process.



⑤ "Wait Bottle Out" is displayed on touch panel when order is completed.





⑥ Put out each cup one by one and fill each bottle with tablets.



⑦ Order status changed to "Checked by". Once the user confirms, then it can attach the applicable label on each bottle.



⑧ Scan the barcode of completed bottle label (order number) for verification by clicking "Checked By" on touch panel. Finally, bottle dispensing for each order is completed.

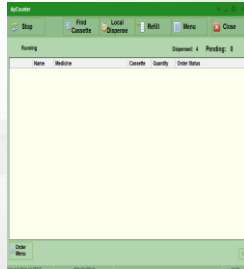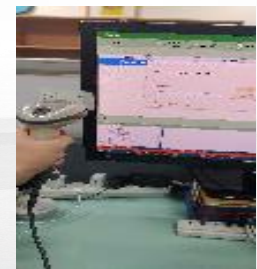

4. CRETEM MTC-30 Product (Max Combination)

MTC-30 1 starter + 2 more extenders



MTC-30 product is semi-vial packager and extendable to max. 90 cassettes.

And one S/W can control all, but, it can't control if the machine is separated.

5-1. MTC-30 Optional Peripheral Device & Specification



1. Touch Panel

- Recommended Product: Microsoft Surface Go 2 LTE or higher
- Standing screen display size: 10.5 inches
- Screen Resolution: 1,920*1,280 Pixels
- Processor: 1.1 GHz
- RAM: 8 GB
- Hard Drive: Minimum 64 GB
- Card Description: Integrated
- Average Battery Life (in hours): 10 hours
- Operating System: Windows 10 Professional
- Platform: Net Framework 3.5
- Communication Device: USB to RS232C

5-2. MTC-30 Optional Peripheral Device & Specification



2. Barcode Scanner

- Recommended Product: ZEBRA SYMBOL-DS2208
- Manufacturer: ZEBRA
- PC Interface Available:
 - 1) PS/2 Method
 - 2) USB Method
 - 3) RS232C Method (Additional Adaptor needed)
- Scan type: 1D/2D, PDF417
- Illuminant: Aiming Pattern – Linear 624nm Amber LED
Illumination – 645nm Super-Red LEDs
- Interface: 39 degrees(H) * 25 degrees (V)
- Reading distance: 0 ~ 38.6 cm

5-3. MTC-30 Optional Peripheral Device & Specification



3. Label Printer

- . Label printer with network connection which can be installed for USB or Driver is recommended
- . Recommended Product: ZEBRA GK420E
- . Manufacturer: ZEBRA
- . Print Method: Direct Sensible Heat
- . Resolution: 203 dpi (Basic-8DOT)
- . Connection Method: USB, RS-232C (Serial), Ethernet
- . Print Speed: Max 152 mm/sec
- . Paper Standard: Width 19 ~ 108 mm
- . Outside Diameter: 127 mm
- . Operating Temperature: 4.4°C ~ 41°C
- . Printing Width: Max 104mm
- . Memory: 8MB flash, 8MB SDRAM
- . Printing Length: Max 991 mm
- . Product Size: 201(L) * 171 (W) * 152 (H) mm
- . Ethernet: LAN Port(RJ45)
- . USB: Universal Serial USB Port
- . Parallel: 36 Pin Parallel Port
- . RS232C(Serial): 9 Pin Direct Port

THANK YOU

